


# How DyeVert™ PLUS Works

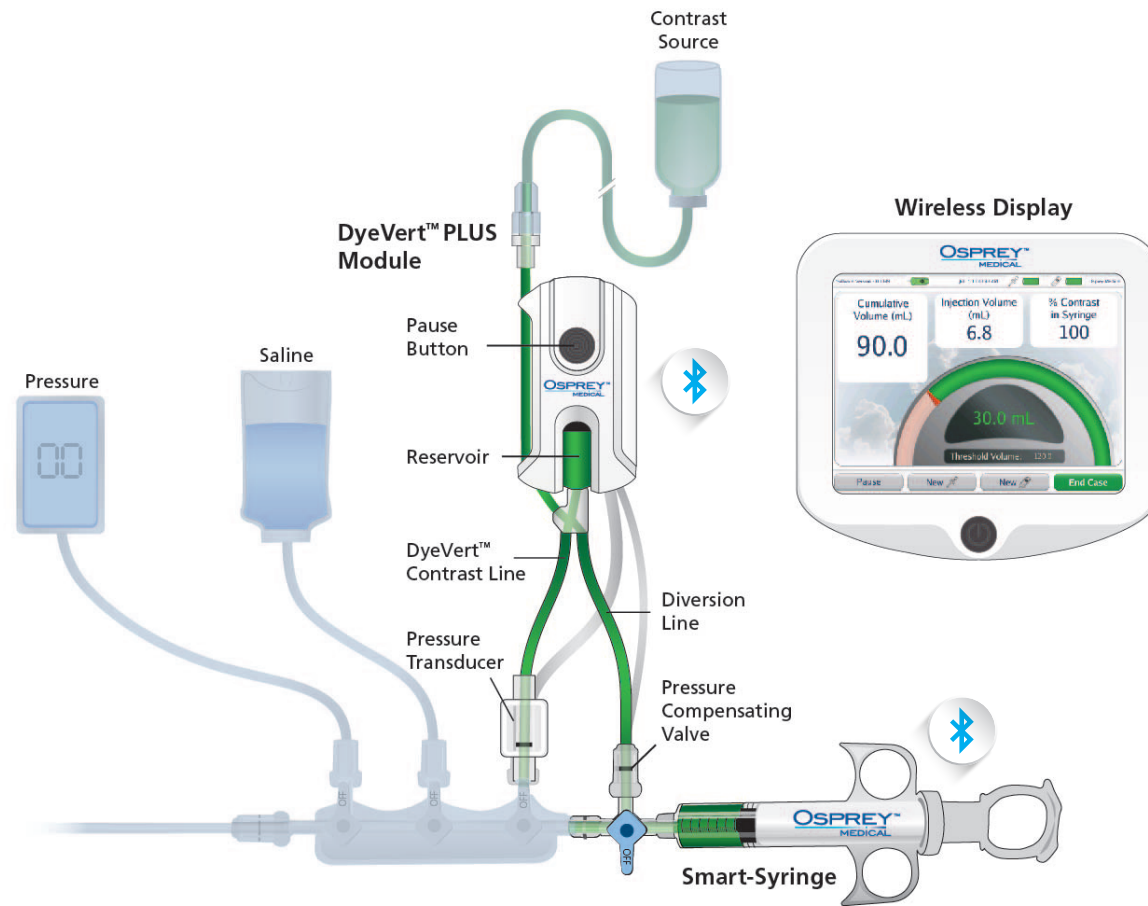

- Interfaces with standard manifold systems
- Easy-to-use, accommodates case flow
- Contrast is diverted on each injection, including puffs
- Reservoir temporarily stores diverted contrast to be used on next injection
- Proprietary Pressure Compensating Valve system automatically self-adjusts for different catheter configurations, contrast types and injection rates

*Indications, contraindications, warnings, and instructions for use can be found in the product labeling supplied with each device. Risks may include air emboli and infection.*

*CAUTION: Federal (USA) law restricts these devices to sale by or on the order of a physician.*

**DYEVERT™**  
Contrast Reduction System **PLUS**